

stryker® INSTRUMENTS

BioPrep Bone Preparation System

Bone Preparation Technique



Femoral Canal Brush-

thoroughly cleans intramedullary canal of loose cancellous bone and other debris.



Cement Restrictor & Inserter-

three ring design provides superior cement fixation and pressurization.




Femoral Sponge- suction feature effectively absorbs blood and fluid from the medullary canal.



Cement Application-

Stryker Dual Speed Gun provides fill and pressurize mode to assure optimal cement-to-bone interface.

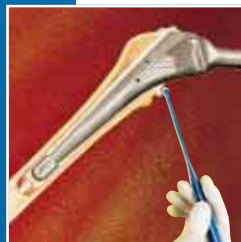

Cement Sculps-

remove surrounding soft cement with no risk of scratching the prosthesis.



206-710

BioPrep Bone Preparation Kit includes:

- Femoral Sponge with Suction Attachment
- Disposable Cement Sculps (2 per kit)
- Femoral Canal Brush
- Universal Restrictor with Disposable Inserter

Convenient

- The entire BioPrep System is available in conveniently prepackaged preparation kit or individually depending on the need and technique that is preferred.

Complete


- The BioPrep System offers three sizes of restrictors as well as a full array of instrumentation for drying and cleansing the femoral canal.

Cost Effective

- Stryker, the proven leader in cement mixing, also offers the BioPrep Total Joint Kit, combining the best products for cement preparation and cement application all in one package. This kit saves time, inventory space and money.

BioPrep

Bone Preparation System

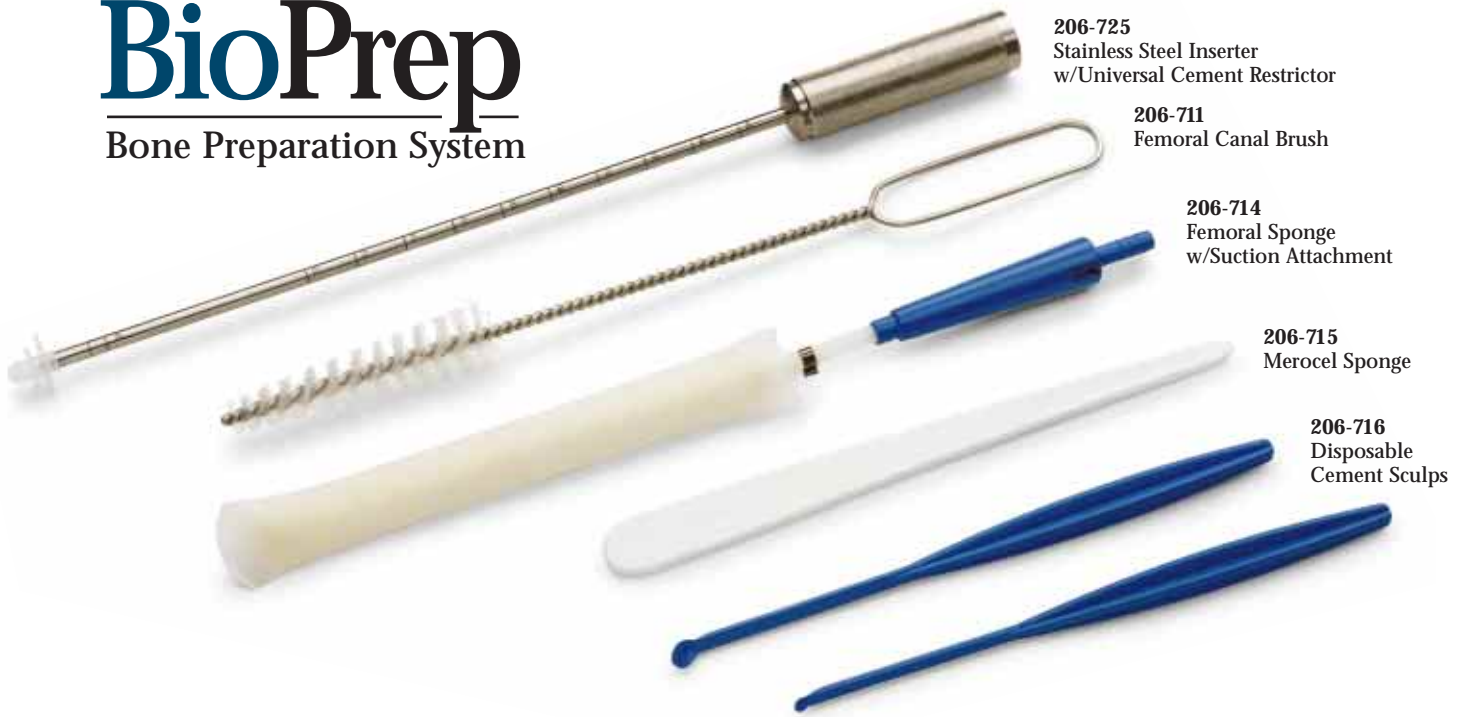

206-725
Stainless Steel Inserter
w/Universal Cement Restrictor

206-711
Femoral Canal Brush

206-714
Femoral Sponge
w/Suction Attachment

206-715
Merocel Sponge

206-716
Disposable
Cement Sculps

Bone Preparation and Cement Application Kits (Sterile, 6 kits per package)

306-705 Advanced Cement Mixing/
BioPrep™ Cement Application Kit
with Medium Pressurizer

- Includes one each of the following:
- Advanced Cement Mixing Bowl w/180 gram Cement Cartridge
 - 200mm Breakaway Femoral Nozzle
 - Universal Cement Restrictor w/ Disposable Inserter
 - Medium Femoral Pressurizer

306-708 Advanced Cement Mixing/
BioPrep™ Total Joint Kit

- Includes one each of the following:
- Advanced Cement Mixing Bowl w/180 gram Cement Cartridge
 - 200mm Breakaway Femoral Nozzle
 - Universal Cement Restrictor w/Disposable Inserter (25mm O.D.)
 - Small Cement Restrictor (18.5mm O.D.)
 - Disposable Cement Sculps (Sterile, 2 per kit)
 - Femoral Canal Brush
 - Femoral Canal Sponge w/Suction Attachment
 - Medium Femoral Pressurizer

206-710 BioPrep™ Bone Preparation Kit

- Includes one each of the following:
- Universal Cement Restrictor w/ Disposable Inserter (25mm O.D.)
 - Small Cement Restrictor (18.5mm O.D.)
 - Disposable Cement Sculps (Sterile, 2 per kit)
 - Femoral Canal Brush
 - Femoral Canal Sponge w/Suction Attachment

Ordering Information

BioPrep™ Bone Preparation System Disposables

(Sterile, 6 per package unless otherwise noted)

- 206-711 Femoral Canal Brush
- 206-714 Femoral Sponge with Suction Attachment
- 206-715 Merocel Femoral Sponge
- 206-716 Disposable Cement Sculps (Sterile, 2 per pouch; 6 pouches per package)
- 206-725 Stainless Steel Inserter (1 per package)

Additional Items

BioPrep™ High Vacuum Cement Mixing Kits

- 206-705 High Vacuum/BioPrep Cement Application Kit
Each Kit Includes One Each of the Following Sterile Items:
 - High Vacuum Mixing Bowl with 180 Gram Cement Cartridge
 - 200mm Breakaway Femoral Nozzle
 - Universal Cement Restrictor with Disposable Inserter
- 206-708 High Vacuum/BioPrep Total Joint Kit
Each Kit Includes One Each of the Following Sterile Items:
 - High Vacuum Mixing Bowl with 180 Gram Cement Cartridge
 - 200mm Breakaway Femoral Nozzle
 - Universal Cement Restrictor with Disposable Inserter
 - Disposable Cement Sculps (Sterile, 2 per kit)
 - Femoral Canal Brush
 - Femoral Sponge with Suction Attachment
 - Medium Femoral Pressurizer

Cement Hardware and Accessories

- 206 Manual High Vacuum Pump
- 206-500 Automatic High Vacuum Foot Pump
- 206-600 Dual Speed Cement Injection Gun
- 206-512 200mm Breakaway Femoral Nozzle
- 206-515 Thin Flexible Nozzle
- 206-522 Acetabular Nozzle
- 206-535 350mm Extra Long Femoral Nozzle
- 206-528 Acetabular Pressurizer
- 206-556 Flat Proximal Pressurizer
- 206-545 Femoral Canal Pressurizer without Hub, Small
- 206-546 Femoral Canal Pressurizer without Hub, Medium
- 206-547 Femoral Canal Pressurizer without Hub, Large

stryker®
INSTRUMENTS

For more information, contact your local Stryker Sales Representative.

4100 East Milham Ave
Kalamazoo, MI 49001
(800)253-3210
(616)323-7700